



# Glove Tester

Device for glove integrity testing.



#### **Dimensions Glove Tester**

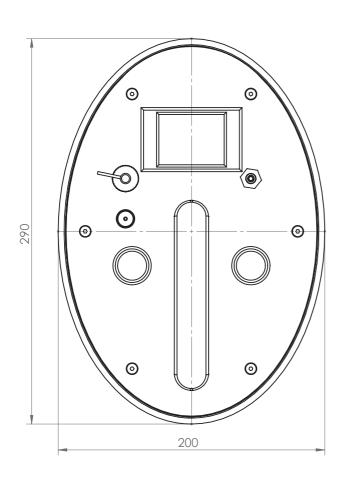

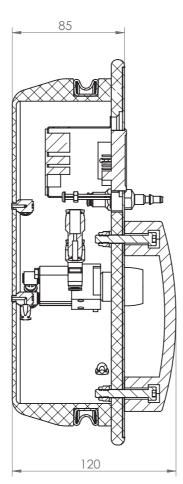

### **Glove Tester**

Dimensions in mm	
Width	200
Height	290
Depth	120

#### Standards

- 01 USB connector
- O Display
- Inflatable seal indicator
- 04 Pressurized air inlet
- **105** Pressure glove inlet/outlet
- 06 Inflatable seal
- 1 Inflatable seal inlet/outlet
- 08 Handle



#### **Operation Readiness**

For operation of the equipment the following must be ensured:

- Pressurised air
- Electric power supply

## Display



#### Description

The tester measures a drop in pressure over time. If the pressure drops over the assigned limit an alarm is displayed. A recording of the displayed values must be performed manually as a part of prepared protocol.